

## CM 23952 - CASEIN ACID HYDROLYSATE (STD) TBL POWDER (Culture Media Ingredients)

### INTENDED USE

Used in antibiotic sensitivity testing media, vaccine preparation media etc.

### PRODUCT SUMMARY AND EXPLANATION

Casein Acid Hydrolysate (Std) TBL is the result of acidic digestion of milk protein by hydrochloric acid. Absence of sulphonamide inhibitors makes Casein Acid Hydrolysate (Std) TBL Powder ideal for preparation of Antibiotic Test Media, like Mueller Hinton Agar and Vaccine Preparation Media as a source of high concentration of free amino acids. It contains all amino acids (except Tryptophan and Cystine which are destroyed during acid hydrolysis) present in milk protein and high sodium chloride content.

### PRINCIPLE

Casein Acid Hydrolysate (Std) TBL is the result of acidic digestion of milk protein-casein by hydrochloric acid. As such, it contains all amino acids (except Tryptophan and Cystine which are destroyed during acid hydrolysis) present in Casein.

### INSTRUCTION FOR USE

Casein Acid Hydrolysate (Std) TBL used as an ingredient in Mueller Hinton Agar.

### QUALITY CONTROL SPECIFICATIONS

Appearance	:	Off white to creamish yellow colour, free flowing powder, having characteristic odour but not pungent smell.
Solubility (2% soln. at 25°C)	:	Soluble in distilled Water, Clear. Insoluble in alcohol.
Clarity (2% Soln. at 121°C)	:	Clear solution. No ppt.
pH (2% Soln. at 25°C)	:	6.5 – 7.5
Loss on drying (at 105°C)	:	NMT – 6.0%
Total Nitrogen (DWB)	:	NLT – 8.5%
α-Amino Nitrogen	:	NLT – 4.5%
Total Ash	:	NMT – 40.0%
Sodium Chloride ( NaCl)	:	NMT – 38.0%
Heavy Metals (Pb)	:	NMT- 20ppm
Indole Test	:	Negative
Microbial test	:	Passes Test

### INTERPRETATION

Cultural Characteristic observed in 1.75% Casein Acid Hydrolysate (Std) TBL used as an ingredient in Mueller Hinton Agar after incubation at 35-37°C for 18-24 hours.

Microorganism	ATCC	Inoculum (CFU/ml)	Growth
Staphylococcus aureus	25923	50-100	Good - Luxuriant
Escherichia coli	25922	50-100	Good - Luxuriant
Pseudomonas aeruginosa	27853	50-100	Good - Luxuriant
Bacillus subtilis	6633	50-100	Good - Luxuriant
Streptococcus pyogenes	19615	50-100	Good - Luxuriant



**PACKAGING:**

Standardpacking is 500gm in HDPE bottle. After packing tightly closed in a dry and well- ventilated place.

**STORAGE**

Store atroom temperature in cool place, Keep HDPE bottle tightly closed in a dry and well-ventilated place. Use before expiry date on label. On opening, product should be properly stored in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use.

Product Deterioration: Do not use product if any contamination, discoloration or other sign of deterioration is found.

**DISPOSAL**


Afteruse,contact a licensed professional waste disposal service to dispose of this material. Dispose of as unused product.

**REFERENCES**

1. J.H. Müller, J. Hinton, Proc. Soc. Exptl. Biol. 48, 330 (1941)
2. G.M. Eliopoulos, et al., Enhancement of cefotaxime and other cephalosporins against Enterococcus faecalis by blood supplemented Mueller-Hinton agar, Diagn. Microbiol. Infect. Dis. 12, 149 (1989);



Quantity

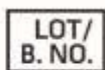

Lot / Batch Number



Temperature Unit



Best Before



QR  
Code



Catalogue No.



Consults Instructions for use :



Manufacturer

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

\*ForLab Use Only

