

CM 20,000 – 2 XYT GROWTH AGAR

INTENDED USE

Optimized formulation for the growth and maintenance of M13 phage or other filamentous ss DNA bacteriophages.

PRODUCT SUMMARY AND EXPLANATION

2XYT Growth Agar is an optimized formulation for the growth and maintenance of M13 phage or other filamentous ss DNA bacteriophages. 2XYT Growth Medium is an optimized formulation for the growth and maintenance of M13 phage or other filamentous ss DNA bacteriophages. This media is 2 times richer than the YT media. This media was originally formulated as a nutritionally enriched growth medium for growth of recombinant strains of Escherichia coli and can also be used for propagation of M13 bacteriophage. It permits larger quantity of phage production without exhausting the host.

COMPOSITION

Ingredients	Gms / Ltr
Tryptone	16.000
Yeast extract	10.000
Sodium chloride	5.000
Agar	15.000

PRINCIPLE

The medium consists of Yeast extract and Tryptone which provide all the required amino acids, nucleotide precursors, vitamins and other metabolites and as a result the cells grow faster in this medium. Sodium chloride provides sodium ions for transport and osmotic balance. This medium contains agar as the solidifying agent.

INSTRUCTION FOR USE

- Dissolve 46.0 grams in 1000 ml purified/distilled water.
- Heat to boiling to dissolve the medium completely.
- Sterilize by autoclaving at 15 psi pressure (121°C) for 15 minutes.
- Mix well and dispense as desired.

QUALITY CONTROL SPECIFICATIONS

- Appearance of Powder : Light yellow coloured, homogeneous, free flowing powder.
- Appearance of prepared medium : Light yellow coloured, clear to slightly opalescent gel forms in Petri plates.
- pH (at 25°C) : 7.0 ± 0.2

INTERPRETATION

Cultural characteristics observed after incubation.

Microorganism	ATCC	Inoculum (CFU/ml)	Growth	Recovery	Incubation Temperature	Incubation Period



Escherichia coli	25922	50-100	Good-luxuriant	>=50%	35-37°C	18-48 Hours
------------------	-------	--------	----------------	-------	---------	-------------

PACKAGING:

Inpacksizeof500 gm bottles.

STORAGE

Dehydrated powder, hygroscopic in nature, store in a dry place, in tightly-sealed containers between 25-30°C and protect from direct sunlight. Under optimal conditions, the medium has a shelf life of 4 years. When the container is opened for the first time, note the time and date on the label space provided on the container. After the desired amount of medium has been taken out replace the cap tightly to protect from hydration.




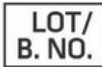








Product Deterioration: Do not use if they show evidence of microbial contamination, discoloration, drying or any other signs of deterioration.

DISPOSAL

Afteruse,prepared plates, specimen/sample containers and other contaminated materials must be sterilized before discarding.

REFERENCES

1. Difcomanual 11th ed., Sparks, MD (1998), 22-23
2. Assubel, F.M., R. Brent, R.E. Kingston, D.D. Moore, J.G. Seidman, J.A. Smith and K. Struhl, Current protocols in molecular biology, vol. 1, Current Protocols, New York, (1994)
3. Davis, L.G., M.D. Dibner and J.F. Battey, Basic methods in molecular biology, Elsevier, new York, (1986).

 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Cataloge Number	 Manufacturer
 Temperature Unit	 EC REP MedMet GmbH Bockstrasse 10 48143 Muenster, Germany Authorized Representative	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only