

## CM 20108 – ANTIBIOTIC ASSAY MEDIUM NO. 19

### INTENDED USE

For microbiological assay of Amphotericin B, Netamycin & Nystatin using *Saccharomyces cerevisiae*.

### PRODUCT SUMMARY AND EXPLANATION

Antibiotic Assay media are used in the performance of antibiotic assays. Grove and Randall have elucidated those antibiotic assays and media in their comprehensive treatise on antibiotic assays. Schmidt and Moyer have reported the use of antibiotic assay medium for the liquid formulation used in the performance of antibiotic assay. This media is prepared according to USP and by FDA. This medium is as per specification of Krishbaum and Arett. Freshly prepared plates should be used for antibiotic assays. Test organisms are inoculated in sterile seed agar precooled to 40-45°C and spread evenly over the surface of solidified base agar.

### COMPOSITION

Ingredients	Gms / Ltr
Peptone	9.400
Yeast extract	4.700
Beef extract	2.400
Dextrose (Glucose)	10.000
Sodium chloride	10.000
Agar	23.500

### PRINCIPLE

Peptone, yeast extract and beef extract provides nutrients and growth factor. Dextrose provides the energy source and sodium chloride maintains the osmotic equilibrium of the medium.

### INSTRUCTION FOR USE

Dissolve 60.0grams in 1000 ml purified / distilled water.

Heat to boiling to dissolve the medium completely.

Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C.

Advice: Recommended in the microbiological assay of Amphotericin B, Candicidin and Nystatin.

### QUALITY CONTROL SPECIFICATIONS

Appearance of Powder	: Cream to yellow homogeneous free flowing powder.
Appearance of prepared medium	: Yellow coloured clear to slightly opalescent gel forms in Petri plates.
pH (at 25°C)	: 6.1±0.2

### INTERPRETATION

Cultural characteristics observed after incubation.

Microorganism	ATCC	Inoculum (CFU/ml)	Growth	Recovery	Antibiotics assayed	Incubation Temperature	Incubation Period



Saccharomyces cerevisiae	2601	50-100	Luxuriant	>=70%	Nystatin	29-31°C	24-48 Hours
Saccharomyces cerevisiae	9763	50-100	Luxuriant	>=70%	Amphotericin B, Candicidin	29-31°C	24-48 Hours

#### PACKAGING:

In pack size of 500 gm bottles.

#### STORAGE

Dehydrated powder, hygroscopic in nature, store in a dry place, in tightly-sealed containers between 25-30°C and protect from direct sunlight. Under optimal conditions, the medium has a shelf life of 4 years. When the container is opened for the first time, note the time and date on the label space provided on the container. After the desired amount of medium has been taken out replace the cap tightly to protect from hydration.


Product Deterioration: Do not use if they show evidence of microbial contamination, discoloration, drying or any other signs of deterioration.

#### DISPOSAL

After use, prepared plates, specimen/sample containers and other contaminated materials must be sterilized before discarding.

#### REFERENCES

1. Grove and Randall, 1955, Assay Methods of Antibiotics Medical Encyclopedia, Inc, New York.
2. Krishbaum A and Areet B, 1967, J. Pharm Sci, 56: 512. 5. Schmidt and Moyer, 1944; J. Bact, 47:199.
3. Tests and Methods of Assay of Antibiotics and Antibiotic containing Drugs, FDA, CFR, 1983. Title 21, part 436, Subpart D, Washington, D.C. U.S Government printing office, paragraphs 436, 100-436, 106 pg 242-259 (April 1).
4. United States Pharmacopoeia 2009. US Pharmacopoeial Convention Inc, Rockville, MD.

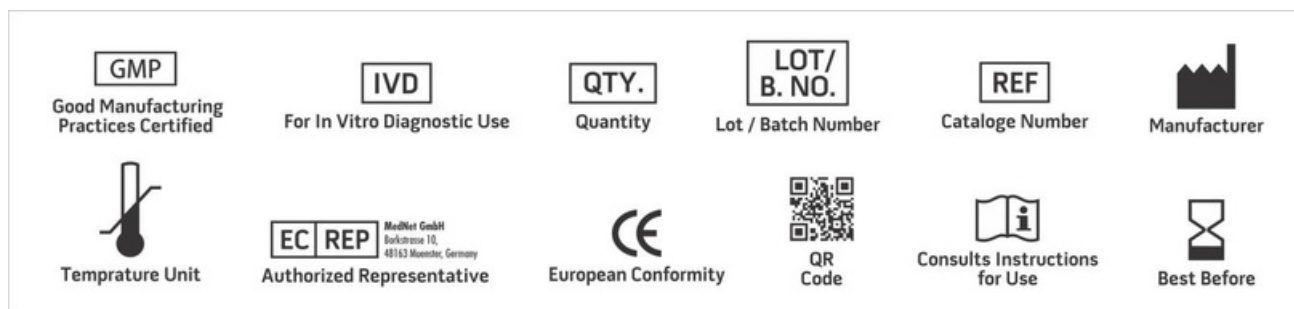

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

\*For Lab Use Only