

CM 20740 – FUNGAL BROTH (MYCOLOGICAL BROTH)

INTENDED USE

For cultivation of fungi.

PRODUCT SUMMARY AND EXPLANATION

Mycological Broth is a basal medium to which antifungal agents may be added for checking their effect on fungi or antibacterials to render them selective for isolation and cultivation of fungi. Mycological Agar is used while working with pathogenic fungi. Earlier media for fungi generally relied on an acidic pH to make the media less suitable for the growth of many bacteria. Fungal Agar is prepared according to the formulation suggested by Huppert and Walker. Mycological Broth is similar in composition to Mycological Agar (Fungal Agar), except agar.

COMPOSITION

Ingredients	Gms / Ltr
Soya peptone	10.000
Dextrose (Glucose)	40.000

PRINCIPLE

The medium consists of Soya peptone in the medium which provides nitrogen, vitamins and minerals necessary to support bacterial growth. Dextrose is a carbon source required for the growth of fungi.

INSTRUCTION FOR USE

Dissolve 50.0 grams in 1000 ml purified/distilled water.
 Heat if necessary to dissolve the medium completely.
 Dispense into tubes or flasks as desired and Sterilize by autoclaving at 15 psi pressure (121°C) for 15 minutes.
 Cool to 45-50°C. For preparing selective media, acidify the media upto pH 3.0-4.0 by the addition of two vials of 10% Lactic Acid Solution .

QUALITY CONTROL SPECIFICATIONS

Appearance of Powder : Cream to yellow homogeneous free flowing powder.
 Appearance of prepared medium : Light amber coloured, clear solution in tubes.
 pH (at 25°C) : 7.0 ± 0.2

INTERPRETATION

Cultural characteristics observed after incubation.

Microorganism	ATCC	Inoculum (CFU/ml)	Growth	Incubation Temperature	Incubation Period
Aspergillus brasiliensis	16404	10-100	Luxuriant	25-30°C	48-72 Hours



Candida albicans	10231	10-100	Luxuriant	25-30°C	48-72 Hours
Lactobacillus acidophilus	11506	50-100	Luxuriant	25-30°C	48-72 Hours
Saccharomyces cerevisiae	9763	10-100	Luxuriant	25-30°C	48-72 Hours
Saccharomyces uvarum	28098	10-100	Luxuriant	25-30°C	48-72 Hours
Staphylococcus aureus subsp. aureus	25923	50-100	Luxuriant	25-30°C	48-72 Hours
Trichophyton mentagrophytes	9533	10-100	Luxuriant	25-30°C	Upto 7 Days

PACKAGING:

Inpacksizeof500 gm bottles.

STORAGE

Dehydrated powder, hygroscopic in nature, store in a dry place, in tightly-sealed containers between 25-30°C and protect from direct sunlight. Under optimal conditions, the medium has a shelf life of 4 years. When the container is opened for the first time, note the time and date on the label space provided on the container. After the desired amount of medium has been taken out replace the cap tightly to protect from hydration.

Product Deterioration: Do not use if they show evidence of microbial contamination, discoloration, drying or any other signs of deterioration.




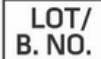








DISPOSAL

After use, prepared plates, specimen/sample containers and other contaminated materials must be sterilized before discarding.

REFERENCES

1. A.J. Clin. Path., 1951, 21: 684.
2. Huppert M., and Walker L. J., 1958, Am. J. Clin. Pathol., 29:291
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
5. Salfinger Y., and Tortorello M.L., 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.



 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Cataloge Number	 Manufacturer
 Temperature Unit	 EC REP Authorized Representative <small>MedNet GmbH Buckstrasse 10 48163 Muenster, Germany</small>	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For LabUse Only

