

CM 23706 – AEROMONAS SELECTIVE SUPPLEMENT

INTENDED USE

For the selective isolation of Aeromonas species.

COMPOSITION

Ingredients	Concentration
Ampicillin	2.500mg

(Per vial sufficient for 500ml medium)

APPLIED IN MEDIUM


AEROMONAS ISOLATION MEDIUM BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml of 50% v/v alcohol. Mix well and aseptically add to 500 ml of sterile, molten, cooled (45-50°C) Aeromonas Isolation Medium Base. Mix well and pour into sterile petri plates.

STORAGE:

Vials should be stored in sealed container at 2-8°C.

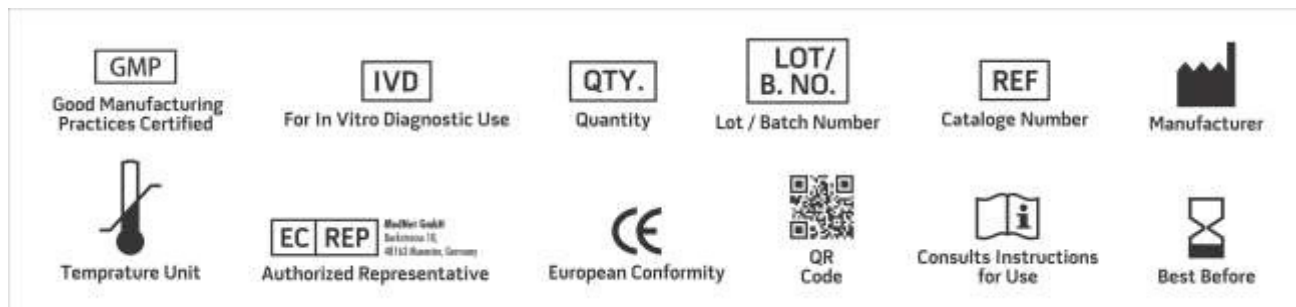

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

