

CM 23709 – AMPICILLIN SUPPLEMENT

INTENDED USE

For isolation of Aeromonas hydrophila from foods.

COMPOSITION

| Ingredients | Concentration |
|-------------|---------------|
| Ampicillin | 10mg |

(Per vial sufficient for 500ml medium)

APPLIED IN MEDIUM

SAAGARBASE (as per APHA)


AEROMONAS STARCH DNA AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml of 50% v/v alcohol. Mix well and aseptically add to 500 ml of sterile, molten, cooled (45-50°C) Aeromonas Isolation Medium Base. Mix well and pour into sterile petri plates.

STORAGE:

Vials should be stored in sealed container at 2-8°C.

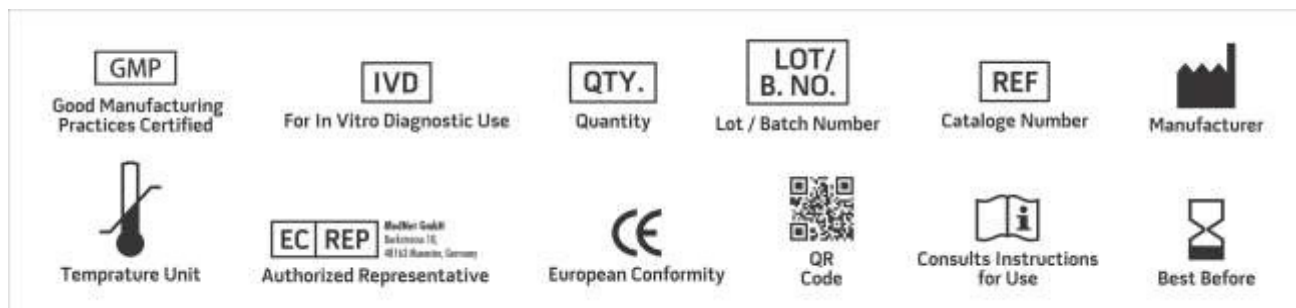

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

