

CM 23713 – BACILLUS SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of Bacillus species.

COMPOSITION

Ingredients	Concentration
Polymyxin sulphate	10mg
Bacitracin	10mg

(Per vial sufficient for 1000 ml medium)

Applied in Medium


CHROMOGENIC BACILLUS AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 10 ml sterile distilled water. Mix well and aseptically add it to 1000 ml sterile, molten CHROMOGENIC BACILLUS AGAR BASE. Mix well and pour into sterile petri plates / tubes.

STORAGE

Vials should be stored in sealed container at 2-8°C.

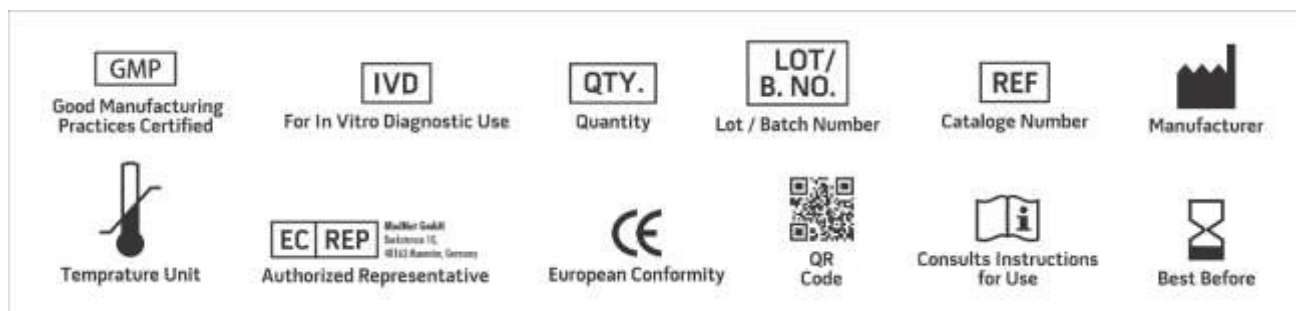

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

