

CM 23719 – BROMO CRESOL PURPLE (15 mg/vl)

INTENDED USE

For selective isolation and enumeration of faecal Streptococci.

COMPOSITION

Ingredients	Concentration
Bromo cresol purple	15mg

(Per vial sufficient for 1000 ml medium)

Applied in Medium


KFSTREPTOCOCCAL AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml of sterile 0.2 N sodium hydroxide and add to 1000 ml of sterile, molten KF STREPTOCOCCAL AGAR BASE, along with 1 vial of 1% TTC Solution. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

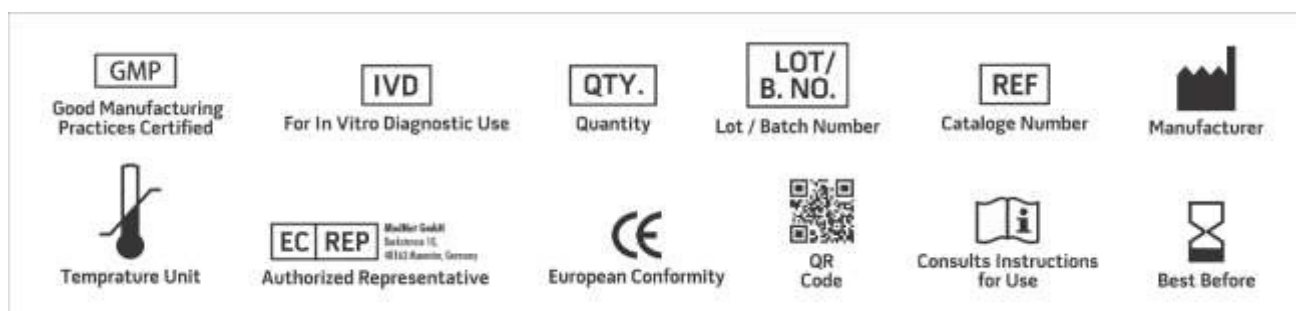

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For LabUse Only

