

CM 23731 – CAMPYLOBACTER SUPPLEMENT-I (Blaser Wang)

INTENDED USE

For selective isolation of Campylobacter species.

COMPOSITION

Ingredients	Concentration
Polymyxin B	1.250IU
Vancomycin	5mg
Trimethoprim	2.500mg
Amphotericin B	1mg
Cephalothin	7.500mg

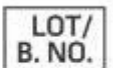
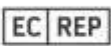
(per vial sufficient for 500ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml sterile distilled water and mix well to dissolve. Avoid frothing of the solution. Aseptically add the contents to 500 ml of sterile, molten, cooled (45-50°C) Blood agar base no. 2 / C. Botulinum isolation agar base / Campylobacter agar base / Columbia blood agar base / Columbia blood agar base (veg.) / Blood agar base no. 2 (ISO 11290-1:2017) / Blood agar base no. 2 with 1.2% agar / Brucella agar base / Brucella agar base (veg.) / Brucella broth base / Campylo thioglycollate medium base (campy-thio medium) / Columbia blood agar base w/ 1% agar with 10% v/v sheep blood or 5-7% v/v laked horse blood and rehydrated contents of 1 vial of Campylobacter Growth Supplement. Mix gently and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

