

CM 23738 – CCDA SELECTIVE SUPPLEMENT

INTENDED USE

For the selective cultivation of Campylobacter species.

COMPOSITION

Ingredients	Concentration
Cefoperazone	16mg
Amphotericin B	5mg

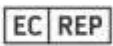
(Per vial sufficient for 500 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, cooled (45-50°C) BLOOD FREE CAMPYLOBACTER BROTH BASE.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative Muller-Guttmann Buckhorn 10, 40143 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*ForLabUseOnly