

CM 23742 – CFC SUPPLEMENT

INTENDED USE

For the selective isolation of *Pseudomonas* species.

COMPOSITION

Ingredients	Concentration
Cetrimide	5mg
Fucidin	5mg
Cephaloridine	25mg

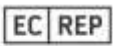
Applied in Medium (Per vial sufficient for 500 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 1 ml 95% Ethanol + 1 ml 0.2N NaOH solution. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) PSEUDOMONAS AGAR BASE. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative MedWet GmbH Buckhorn 16 40162 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only