

CM 23755 – CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT

INTENDED USE

For the isolation, detection, and enumeration of *Listeria monocytogenes*.

COMPOSITION

Ingredients	Concentration
Ceftazidime	2mg
Amphotericin B	2.500mg

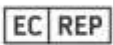
(Per vial sufficient for 500ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 5 ml of sterile distilled water. Mix well and aseptically add to 500 ml sterile, molten, cooled (45-50°C) CHROMOGENIC L-MONO LISTERIA DIFFERENTIAL AGAR / CHROMOGENIC LISTERIA AGAR BASE (ISO 11290-1 & 2:2004, ISO 11133:2014) / CHROMOGENIC LISTERIA AGAR BASE (Modified). Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative MedWit GmbH Buckhorn 10, 40143 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only