

CM 23762 – CLOSTRIDIUM PERFRINGENS SUPPLEMENT

INTENDED USE

For selective isolation of Clostridium perfringens.

COMPOSITION

Ingredients	Concentration
D-Cycloserine	200mg
4-Methylumbelliferyl-phosphate disodium salt	50mg

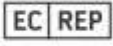
(Per vial sufficient for 500 ml medium)

INSTRUCTION FOR USE

Rehydrate the content of 1 vial aseptically with 3 ml sterile distilled water. Mix well and aseptically add to 500 ml of sterile, molten, cooled (45-50°C) PERFRINGENS AGAR BASE (TRYPTOSE SULPHITE CYCLOSERINE AGAR BASE) (ISO 7937: 2004, ISO 14189:2013) / TRYPTONE YEAST SODIUM SULPHITE AGAR BASE (ISO 14189:2013). Mix well and pour into the sterile petri plates.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 EC REP Authorized Representative MedWit GmbH Buckhorn 16 48163 Münster, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For LabUse Only