

CM 23770 – CHROMOGENIC EC O157:H7 SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation and detection of Escherichia coli O157:H7 from food samples.

COMPOSITION

Ingredients	Concentration
Novobiocin	10mg
Potassium tellurite	1mg

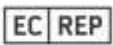
(Per vial sufficient for 1000ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 10 ml of sterile distilled water. Mix well and add aseptically to 990 ml sterile, molten, cooled (45-50°C) CHROMOGENIC EC O157:H7 AGAR (ISO 16654:2001) / CHROMOGENIC ENRICHMENT BROTH BASE FOR EC O157:H7 / BUFFERED PEPTONE WATER. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 EC REP Authorized Representative Machner GmbH Buckhorn 10, 40142 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only