

CM 23776 – ENTEROCUCUS SELECTIVE SUPPLEMENT

INTENDED USE

Enumeration of intestinal Enterococci from surface and wastewater.

COMPOSITION

Ingredients	Concentration
Nalidixic acid	25mg
Thallous acetate	200mg

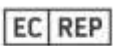
(per vial sufficient for 100ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 5 ml of sterile distilled water. Mix well and aseptically add it to sterile, cooled (45-50°C) 100 ml of MUG SF BROTH BASE / MUD SF BROTH BASE along with 1 ml of sterile 1% T.T.C. Solution. Mix well and dispense as desired.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 EC REP Authorized Representative Müller-Guttmann Biochemie AG 48143 Münster, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.