

CM 23780 – FRASER SUPPLEMENT

INTENDED USE

For selective isolation and enumeration of *Listeria monocytogenes* from food, animal feeds, etc.

COMPOSITION

Ingredients	Concentration
Ferric ammonium citrate	250mg

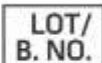
(per vial sufficient for 500 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 1-2 ml of sterile distilled water. Mix well and for primary enrichment aseptically add 1 vial of 500 ml, sterile, cooled (45-50°C) FRASER BROTH BASE, MODIFIED (HALF FRASER BROTH) / FRASER BROTH BASE (ISO 11290-1: 2017) / HALF FRASER BROTH BASE (ISO 11290-1:2017) / FRASER BROTH BASE along with rehydrated contents of 1 vial of Fraser Selective Supplement. Mix well and dispense as desired.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative WuXi AppTec GmbH Helmholtzstr. 10 40113 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.