

CM 23786 – GENTA-OXY SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of vancomycin-resistant enterococci.

COMPOSITION

Ingredients	Concentration
Oxytetracycline	50mg
Gentamicin	20mg

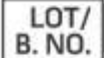
(per vial sufficient for 500 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) OXYTETRA GLUCOSE YEAST AGAR BASE (OGYE AGAR BASE) ISO 6611:2004. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only