

CM 23796 – KL VIRULENCE ENRICHMENT (20 ml/vl)

INTENDED USE

For cultivation and in-vitro toxicity testing of *Corynebacterium diphtheriae*.

COMPOSITION

Ingredients	Concentration
Casein acid hydrolysate	10g
Glycerol	10ml
Polysorbate 80	10ml

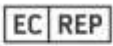
(per vial sufficient for 100 ml medium)

INSTRUCTION FOR USE

Warm up the refrigerated contents of 1 vial to 50°C and aseptically add 2 ml in 100 mm sterile petri plate along with 0.5 ml of 1% Potassium Tellurite Solution. Quickly add 10 ml sterile molten, cooled (45-50°C) DIPHTHERIA VIRULENCE AGAR BASE into the plate. Mix well.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative MedWet GmbH Buckhorn 10, 48147 Münster, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only