

CM 23805 – LACHICA SUPPLEMENT

INTENDED USE

For selective isolation of *Aeromonas hydrophilla*.

COMPOSITION

Ingredients	Concentration
Ampicillin	10mg

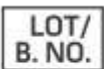
(Per vial sufficient for 1000 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 2 ml of 50% v/v alcohol. Mix well and aseptically add to 1000 ml sterile, molten LACHICA'S MEDIUM BASE. Mix well and pour into sterile petri plates.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

