

CM 23846 – MUELLER TELLURITE SERUM (25 ml/vl)

INTENDED USE

For isolation, cultivation and differentiation of *Corynebacterium diphtheriae*.

COMPOSITION

Ingredients	Concentration
Potassium tellurite	100mg
Calcium D pantothenate	0.050mg
Horse serum	12.500ml
Sodium lactate solution	10ml
Ethanol	2.500ml

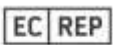
(Per vial sufficient for 1000 ml medium)

INSTRUCTION FOR USE

Warm up the refrigerated solution to room temperature. Add the contents of 1 vial to 975 ml of sterile, molten, cooled (45-50°C) MUELLER TELLURITE AGAR BASE. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative Mueller GmbH Buckhorn 10, 40147 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only