

CM 23854 – NEO ENRICHMENT SELECTIVE SUPPLEMENT

INTENDED USE

Recommended for the selective isolation of *Listeria* species from food samples in 24 hours.

COMPOSITION

Ingredients	Concentration
Ceftazidime	15mg
Nalidixic acid	10mg

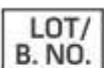
(per vial sufficient for 500ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 5 ml of sterile distilled water and aseptically add to 500 ml of sterile, cooled Neo Enrichment Broth Base. Mix well and dispense in sterile test tubes.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

