

CM 23783 – POLYMYXIN B SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of various microorganisms.

COMPOSITION

Ingredients	Concentration
Polymyxin B sulphate	50000 Unit

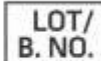
(per vial sufficient for 500/1000 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 2 ml sterile distilled water. Mix well and aseptically add it to 475 ml of sterile, molten BACILLUS CEREUS AGAR BASE or to 450 ml of KG AGAR BASE / BACILLUS CEREUS SELECTIVE AGAR BASE (MYP) (ISO 7932:2004) / MYP AGAR BASE (PHENOL RED EGG YOLK POLYMYXIN AGAR BASE) (IS : 5887 (Part VI) 1976, reaffirmed 2005) / MYP AGAR BASE (PHENOL RED EGG YOLK POLYMYXIN AGAR BASE) (VEG.) / MANNITOL EGG YOLK POLYMYXIN AGAR (ISO 7932:2004, ISO 21871:2006) along with 25 ml/50 ml Egg Yolk Emulsion to make a total volume of 500 ml or to 500 ml of SDS AGAR (SODIUM DODECYL SULPHATE POLYMYXIN SUCROSE AGAR) (as per APHA) / SALT POLYMYXIN BROTH BASE / SALT POLYMYXIN BROTH BASE (ISO 8914:1990) / SOYBEAN CASEIN DIGEST MEDIUM BASE W/O POLYMYXIN / CHROMOGENIC BACILLUS AGAR. Mix well and pour into sterile petri plates / tubes.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative Mallikarjun Gokul Bachchan 16, 10113 Karamba, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only