

CM 23877 – POTASSIUM TELLURITE 1% (1 ml /v)

INTENDED USE

For the selective isolation of Staphylococci and Corynebacteria.

COMPOSITION

Ingredients	Concentration
Potassium tellurite Concentrate	1.100ml

(To achieve 1% solution dilute the contents in 8.9 ml sterile distilled water)


INSTRUCTION FOR USE

Warm up the refrigerated contents of one vial to room temperature. Add aseptically 8.9 ml sterile distilled water, mix well and add in sterile, molten, cooled (45-50°C) AGAR MEDIUM O (BAIRD PARKER AGAR BASE) (as per USP/EP/BP) / BAIRD PARKER AGAR BASE (as per EP/IP/BP) / BAIRD STAPHYLOCOCCUS ENRICHMENT BROTH BASE / VOGEL-JOHNSON AGAR BASE W/O TELLURITE (V. J. AGAR BASE) / VOGEL JOHNSON AGAR BASE W/O TELLURITE (V. J. AGAR BASE) (VEG.) / VOGEL JOHNSON AGAR MEDIUM (as per USP) / VOGEL JOHNSON AGAR MEDIUM (as per IP) / MYCOPLASMA BROTH BASE W/ CV (PPLO BROTH BASE W/ CV) / TPEY AGAR BASE / TELLURITE GLYCINE AGAR BASE / CHOLERA MEDIUM BASE / TM 707



STORAGE

Vialsshouldbe stored in sealed container at 2-8°C.

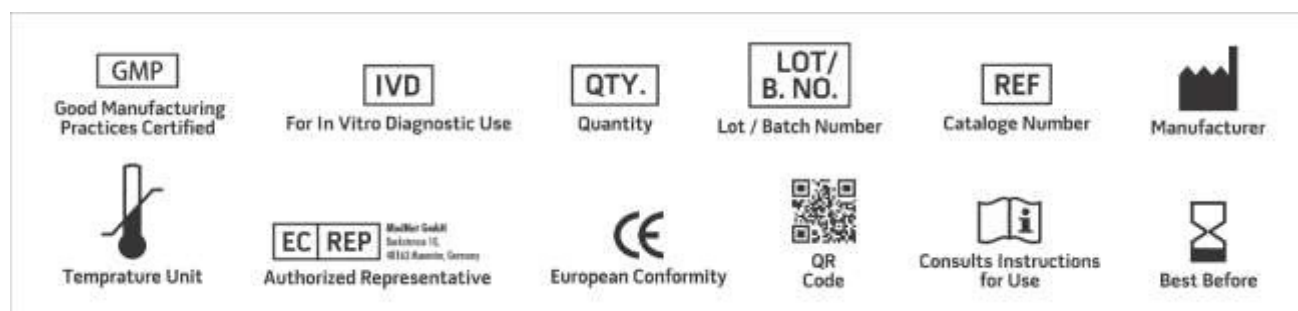

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For LabUse Only