

CM 23886 – SALENRICH SELECTIVE SUPPLEMENT

INTENDED USE

For two-step enrichment of sub-lethally injured Salmonellae from food and feeds.

COMPOSITION

Ingredients	Concentration
Potassium tetrathionate	5g
Ox bile	2g
Brilliant green	0.0175g

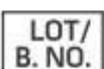
(per vial sufficient for 250 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 25 ml sterile distilled water. Aseptically add 1 ml each of reconstituted SalEnrich Selective Supplement to 10 ml of SALENRICH BROTH at the secondary selective enrichment stage. Allow to stand for 30 minutes. Shake vigorously and then incubate further for 18-22 hours at 35-37°C. For further detection of Salmonella, streak a sample onto appropriate selective media.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only