

CM 23887 – SALMONELLA SELECTIVE ENRICHMENT SUPPLEMENT

INTENDED USE

Recommended for selective isolation and differentiation of Salmonella species from coliforms by chromogenic method.

COMPOSITION

Ingredients	Concentration
Novobiocin	6mg

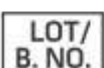
(per vial sufficient for 500ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 5ml of sterile distilled water. Mix gently to dissolve the contents completely. Aseptically add the rehydrated contents of one vial to 500ml of sterile, cooled (45-50°C) of SALMONELLA SELECTIVE ENRICHMENT BROTH BASE. Mix well and dispense as desired.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

