

CM 23890 – SHIGELLA SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of Shigella species.

COMPOSITION

Ingredients	Concentration
Novobiocin	0.550mg

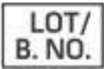
(per vial sufficient for 1000ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 5 ml sterile distilled water and aseptically add it to 1000 ml of sterile, molten, cooled (45-50°C) SHIGELLA BROTH BASE. Mix well and dispense as desired.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only