

CM 23895 – STERILE CHARCOAL SUPPLEMENT FOR LEGIONELLA AGAR

INTENDED USE

Recommended for growth of Legionella species.

COMPOSITION

Ingredients	Concentration
Activated charcoal	2 gm

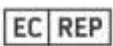
(Per vial sufficient for 500 ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 10 ml of sterile distilled water to dispense the contents evenly. Aseptically add to 430 ml sterile, molten, cooled (45-50C) Legionella Agar Base along with 1 vial containing 50 ml of BCYE. Growth supplement. Aseptically add 10 ml of sterile distilled water to bring volume to 500 ml when no selective supplement is added. Mix well and pour into sterile Petri plates. The final pH of the medium will be 6.9 ± 0.2 .

STORAGE

Vial should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative MedWet GmbH Seckestraße 15 40162 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

