

CM 23896 – STREPTO SUPPLEMENT

INTENDED USE

For selective cultivation Streptococcs species.

COMPOSITION

Ingredients	Concentration
Nalidixic acid	3.750mg
Neomycin sulphate	1.060mg
Polymyxin B sulphate	8500Units

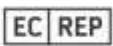
(per vial sufficient for 500ml medium)

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2-5 ml of sterile distilled water. Aseptically add the contents to 500 ml of sterile, molten, cooled (45-50°C) COLUMBIA BLOOD AGAR BASE / COLUMBIA BLOOD AGAR BASE (VEG.) / TM 2044 - COLUMBIA BLOOD AGAR BASE W/ 1% AGAR / BLOOD AGAR BASE No. 2 / BLOOD AGAR BASE NO. 2 (ISO 11290-1:2017) / BLOOD AGAR BASE NO. 2 WITH 1.2% AGAR. The medium may be enriched with 5-7% v/v sterile defibrinated blood. Mix gently and pour into sterile Petri plates.

STORAGE

Vial should be stored in sealed container at 2-8°C.

 Good Manufacturing Practices Certified	 For In Vitro Diagnostic Use	 Quantity	 Lot / Batch Number	 Catalogue Number	 Manufacturer
 Temperature Unit	 Authorized Representative MedWet GmbH Buckhorn 10 48143 Münster, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only

