

CM 23904 – TRICHODERMA HARZIANUM SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation of *Trichoderma harzianum*.

COMPOSITION

Ingredients	Concentration
Chloramphenicol	0.250g
Streptomycin	0.090g
Quintozene	0.200g
Propamocarb	0.926g
Distilled water	40ml

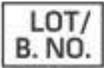
(per vial sufficient for 500ml medium)

INSTRUCTION FOR USE

Warm up the contents of the refrigerated vial to room temperature. Aseptically add to 960ml of sterile, cooled TRICHODERMA HARZIANUM SELECTIVE AGAR BASE. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 Best Before	 QTY. Quantity	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 LOT/ B. NO. Lot / Batch Number	 Consults Instructions for Use	 QR Code	

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only